

## ACCESSORIES

### Optics

Incident beam mask, Soller slits, divergence slit, direct beam absorber, set of alignment slits

### Sample holders

Standard holders (deep and shallow cavities, perforated), Standard Reference material holder, glass holders

### K-beta filter

Ni 0.02 mm filter for diffracted beam

## APPLICATIONS

- Science and education
- Geology and mining
- Cement industry
- Chemical industry
- Oil-gas exploration
- Forensic expertise
- Automotive industry
- Pharmaceutical and cosmetics industries
- And others... (archeology, gemology, food industry etc.)




Start Your Journey with the  
**TELLUS**  
XRD series

### Standard Reference material

NIST 1976c

### Optional Accessories

- Anti-scattering slit
- Zero diffraction plate
- Rotating sample holder
- Sample changer on 8 pcs
- Knife edge - collimator
- Heating chamber
- Variable divergence slit (VDS)
- Holder for the bulk samples

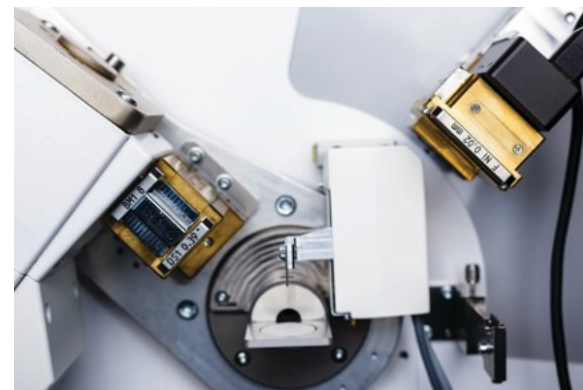

## TECHNICAL DATA SHEET

Generator power	600/300 W (on request)
Target	Cu (Cr, Fe, Co, Mo - optional)
Cooling system	Internal water cooling: closed-circuit
Goniometer	$\theta$ - $\theta$ vertical, radius 150 mm
Scanning range	$-6^\circ \div +154^\circ$ ( $2\theta$ )
Minimal step	$0.001^\circ$ ( $2\theta$ )
Positioning accuracy	$\pm 0.01^\circ$ ( $2\theta$ )
Scanning speed	0.01-600°/min
Dimensions (W x D x H)/Weight	700x700x820 mm/115 kg
Personal computer	PC with OS Windows
Interface	USB/Ethernet

**LINEV SYSTEMS**  
United for Innovation

linevsystems.com  
info@linevsystems.com

© LINEV Systems. All rights reserved.  
This material is provided for informational purposes only and does not modify or replace the applicable terms and conditions. In the event of any conflict, the terms and conditions shall govern and control. Product design and specifications may change without notice. Refer to the official technical proposal for current and binding specifications. Unauthorized use, reproduction, disclosure, or export is prohibited.

**LL-C**  
Certification  
9001:2015

0912-AA12032026



**TELLUS**  
XRD series

Multi-function solution with exceptional materials research & analysis speed, reliability and reproducibility X-ray diffractometer that meets the highest requirements for this class of instruments and even exceeds expectations

## VERY USEFUL AND POWERFUL TECHNIQUE FOR

RESEARCH AND EDUCATION

QUALITY CONTROL

ENVIRONMENTAL, GEOLOGY AND MINERALOGY

PHARMACEUTICALS, DRUG DISCOVERY

CEMENT AND OIL INDUSTRY

FORENSIC SCIENCE

CHEMICAL PROCESS DEVELOPMENT, ANALYTICAL CHEMISTRY

## X-RAY POWDER DIFFRACTOMETER

DESKTOP ANALYTICAL INSTRUMENTS

**LINEV SYSTEMS**  
United for Innovation



From Classroom to Advanced Scientific Research:  
**Precision You Can Trust**



## SOLUTION WITH OPPORTUNITIES OF ANALYTICAL TOOL AT A LOW COST

TELLUS is a compact and advanced tool designed for X-ray diffraction analysis. It is used for researching and identifying new compounds, as well as for quality control and determining the composition of different crystalline substances such as solid samples and powders.

TELLUS is a valuable asset in academic laboratories, enabling researchers to conduct precise studies on synthesized samples. It also aids students in quickly and effectively learning X-ray diffraction techniques.

The diffractometer can be equipped with either 1D or 2D photon counting Si detector with a pixel width of 50  $\mu\text{m}$  providing fast and high-resolution measurements. Fast, robust and maintenance free detectors with 100% efficiency for X-ray radiation with minimum loss of intensity and maximum linear range. The high precision class goniometer together with small pixel width of detector ensure the angle accuracy of better than  $\pm 0.01^\circ$  (2 $\theta$ ). The instrument with proper optic components provides the peak width at least 0.05 $^\circ$  2 $\theta$  over the full angular range.

The software interface with predefined experimental scenarios makes the measurement process clear even for new users, though the experts in XRD can create their own protocol. The software for analysis is based on the advanced analytical methods to make the quantitative and qualitative analysis of samples more precise.

Both the detector and the X-ray source have a long lifetime, what significantly reduces the maintenance costs and saves the budget.

	Mythen2 R 1D	MiniPIX TPX 3
<b>Pixel size</b>	50 $\mu\text{m}$ by 8 mm	55 x 55 $\mu\text{m}$
<b>Type / number of pixel</b>	Linear 1D / 640	Matrix 2D / 250x250
<b>Energy range (keV)</b>	4 - 40	2 - 60
<b>Detection mode</b>	Single threshold	Two-threshold
<b>Frame rate (Hz)</b>	100	16
<b>Cooling</b>	Air	Air
<b>Module weight (g)</b>	100	41

## KEY FEATURES

- Compact design
- Theta/theta geometry
- High effective detector
- Fast measurements
- High signal/noise ratio
- Long lifetime X-ray tube and detector
- No external cooling
- Safety guaranteed
- Analytical software



## SOFTWARE PACKAGE

### TELLUSCon

for system control and data acquisition

- Interface optimized for advanced and non-experienced users
- Measurement by steps or continuous scan
- Predefined measurement scenarios
- Data preprocessing
- Measurement report

### SPECIALISED SOFTWARE

for advanced data analysis

- Data preprocessing (virtual monochromator, Lorentz-polarization factors correction, absorption correction, smoothing, background definition)
- Peak search and its identification
- Qualitative phase analysis
- Quantitative analysis with RIR and internal standard calibration methods
- Refinement of lattice parameters by whole pattern-fitting method (Rietveld, Pawley, Le Bail)
- Analysis report

### Databases

- COD (open-access collection)
- ICDD PDF-2 (optional)

